# SALTER

## BATHROOM SCALES

WITH (WALL-MOUNTED) MIRROR UNIT

### Instructions and Guarantee

---

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---

**PERSONENWAAGE**

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---

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---

**BILANCIA DA BAGNO**

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---

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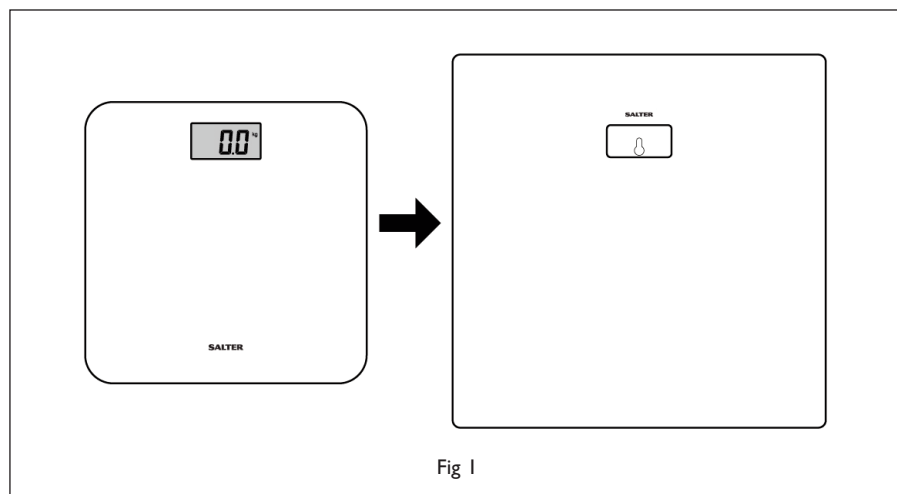
---

Your quality Salter Mirror Scale has been designed especially for your convenience. The wall mounted mirror unit provides a practical dual-purpose storage solution for your scale. When not in use the scale becomes a useful clock.

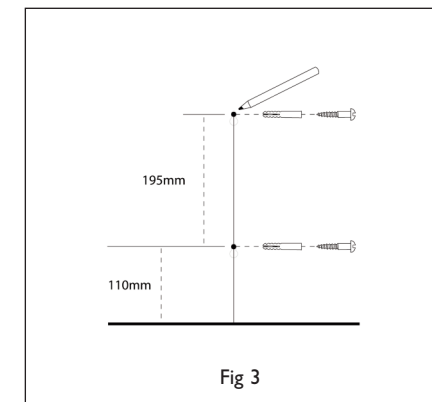
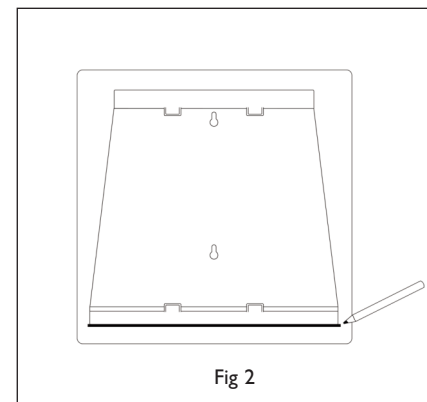
### FIXING YOUR MIRROR UNIT TO THE WALL

**Note:** To fix the mirror unit to your bathroom wall you will need a drill, a screwdriver, a pencil and a spirit level.

**CAUTION:** Please ensure the wall you wish to fit your mirror unit to does not have any cables, pipes (or similar) behind it.

1. Observing the above caution, choose a suitable location for the mirror unit ensuring the scale can be comfortably accessed. Allow sufficient space around the mirror unit for removing and replacing the scale (Fig 1). The scale may be removed or replaced from either side.

2. With a spirit level, carefully position the mirror unit on the wall at your desired height and mark a horizontal pencil line on the wall from corner to corner along the bottom of the mirror unit (Fig 2).
  3. With a spirit level mark a vertical centreline starting from the middle of the horizontal line you have just drawn (Fig 3).
  4. Mark the screw hole positions at 110mm up and 195mm apart (Fig 3).
  5. Using a suitable drill, drill holes for the wall plugs (The wall plugs provided are suitable for solid walls and require a 6 mm diameter hole. For hollow or other wall types – use a suitable alternative fixing.).
  6. Using a suitable screwdriver fit the screws, but do not fully tighten.
  7. Temporarily position the mirror unit and adjust the screws until the mirror unit is held firmly against the wall (You will need to remove the mirror unit to adjust the screws).
- Tip:** Look through the clock window to locate the screw heads.
8. Fit the mirror unit onto the screws and pull down to ensure it is correctly located the screws sitting in the upper part of the keyholes.
  9. Before storing the scale for the first time check the mirror unit is secure and level.
  10. To remove the scale, reach behind the mirror and locate the vertical bar on the back of the scale glass platform. Slide the scale to the edge of the mirror unit until you can firmly grasp the scale platform. Place your other hand under the bottom edge of the scale platform and carefully remove. To replace the scale simply slide the scale behind the mirror. The scale is guided into position and should rest with the clock clearly visible through the mirror window.

**CAUTION:** For safety, the scale should only be removed from the mirror unit with clean dry hands.

## PREPARING YOUR SCALE

1. Open the battery compartment on the scale underside.
2. Insert batteries (2 x AAA) observing the polarity signs (+ and -) inside the battery compartment.
3. Close the battery compartment.
4. Scale enters clock set mode.

## TO SET CLOCK

1. The scale features a 24 hour clock. To enter clock set mode insert new batteries or press and hold **'SET'** for 3 seconds. The hour digits will begin to flash.
2. Press **'UNIT/UP'** button to set hours. Press and hold to advance quickly.
3. Press **'SET'** to confirm. The minute digits will begin to flash.
4. Press **'UNIT/UP'** button to set minutes. Press and hold to advance quickly.
5. Press **'SET'** to confirm.

## OPERATING YOUR SCALE

1. Select kg, st or lb weight mode by pressing the **'UNIT/UP'** switch on the underside of the scale.
2. Position scale on firm flat surface. (For use on carpet remove the non-slip pads from the scales feet and attach the enclosed carpet feet).
3. Tap the platform centre (vibration activates your scale) and remove your foot.
4. Wait until zero is displayed.
5. Step on and stand very still while the scale computes your weight.
6. Your weight is displayed.
7. Step off.
8. Your weight is displayed for a few seconds then the scale will display the clock.

## WARNING INDICATORS

- '----'** Weight exceeds maximum capacity.
- 'Lo'** Replace battery.

## ADVICE FOR USE AND CARE

Always weigh yourself on the same scale placed on the same floor surface. Do not compare weight readings from one scale to another as some differences will exist due to manufacturing tolerances.

- Placing your scale on a hard, even floor will ensure the greatest accuracy and repeatability.
- Weigh yourself at the same time each day, before meals and without footwear. First thing in the morning is a good time.
- Your scale rounds up or down to the nearest increment. If you weigh yourself twice and get two different readings, your weight lies between the two.
- Clean your scale with a damp cloth. Do not use chemical cleaning agents.
- Do not allow your scale to become saturated with water as this can damage the electronics.
- Treat your scale with care - it is a precision instrument. Do not drop it or jump on it.

## TROUBLESHOOTING GUIDE

If you experience any difficulty in the use of your scale:

- Check the battery is correctly fitted.
- Check that you have selected your preference of stones/pounds, kilograms or pounds.
- Check that the scale is on a flat, level floor and not touching against a wall.
- If when you use the scale there is no display at all or if **'Lo'** is displayed - try a new battery.
- If **'----'** is displayed the scale has been overloaded.

## GUARANTEE

This product is intended for domestic use only. Salter will repair or replace the product, or any part of this product, (excluding batteries) free of charge if within 15 years of the date of purchase, it can be shown to have failed through defective workmanship or materials. This guarantee covers working parts that affect the function of the scale. It does not cover cosmetic deterioration caused by fair wear and tear or damage caused by accident or misuse. Opening or taking apart the scale or its components will void the guarantee. Claims under guarantee must be supported by proof of purchase and be returned carriage paid to Salter (or local Salter appointed agent if outside the UK). Care should be taken in packing the scale so that it is not damaged while in transit. This undertaking is in addition to a consumer's statutory rights and does not affect those rights in any way. For UK Sales and Service contact HoMedics Group Ltd, PO Box 460, Tonbridge, Kent, TN9 9EW, UK. Helpline Tel No: (01732) 360783. Outside the UK contact your local Salter appointed agent.

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