HOMEDICS®

Eye and Temple Massager

IM-100-EU

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

IMPORTANT SAFEGUARDS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- DO NOT place or drop into water or any other liquid.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has been damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, send to a HoMedics Service Centre. All servicing of this appliance must be performed by authorised HoMedics service personnel only.

SAFETY CAUTIONS:

PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

- If you have any concerns regarding your health, consult a doctor before using this appliance.
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time of 15 minutes.

This product contains small parts that may present a choking hazard for young children.

WARNING:

1. Do not use this unit if you have had an eye operation recently. If you suffer from any serious eye disease, consult your doctor before using this unit.
2. Do not wear glasses or contact lenses while using this unit.
3. The eye mask must be worn when turning the power on. This is to keep the air pocket from being damaged.
4. Stop the unit immediately if you feel any physical discomfort or if you detect anything unusual.
5. Eyes must be closed while the unit is in operation.

6. Ensure the eye mask and elastic strap are comfortable.
7. Avoid excessive bending or pulling on the connection between the cable, remote control and eye mask. Doing so will stop the air flow and effect the feel of the massage.

BATTERY CAUTIONS:

- Use only the size and type of batteries specified.
- When installing batteries, observe proper +/- polarities. Incorrect installation of battery may cause damage to the unit.
- Do not mix different types of batteries together (e.g., alkaline with carbon-zinc or old batteries with new ones).
- If the unit is not going to be used for a long period of time, remove batteries to prevent damage due to possible battery leakage.
- Do not dispose of batteries in fire. Batteries may explode or leak.

MAINTENANCE:

Cleaning

Turn off the eye mask before cleaning.

- Wipe clean with a soft, damp cloth.
- Inspect the eye mask, controller and power cable to ensure there is no damage.


Storage

Store in a cool, dry place, away from heat, humidity and pollution.

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

PRODUCT FEATURES:

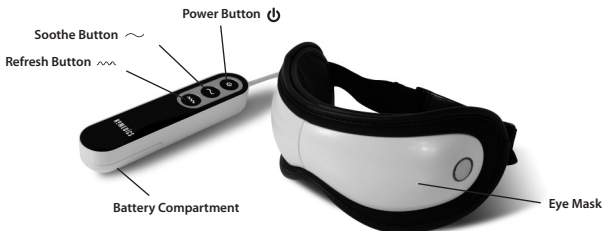INSTRUCTIONS FOR USE:

Inserting the Batteries

1. Remove the battery compartment cover on the controller.
2. Insert 3 'AA' batteries (included) according to the polarity indicators inside the compartment.
3. Replace the battery compartment cover.

Using the Eye Mask

1. Place the eye mask over your eyes.
2. Secure the straps around your head with the clip closure.
3. Press **⏻** power button to turn the eye mask on.
4. The unit will start in soothe mode, the lightest air pressure setting.
5. Press soothe **~** button squeeze-and-hold inflation provides maximum pressure relief for sinuses and temples.
6. Press refresh **~ ~ ~** button gradual full inflation delivers wave-like pressure and relief.
7. Press **⏻** power button to turn the eye mask off.

Note: The eye mask shuts off automatically after 15 minutes. When the eye mask is turned off and back on, it resumes the program last used

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