# HOMEDICS

THERA : P

## Massaging Hot Water Bottle

**Instruction Manual**

---

WTR-BOT-0EU

**2** year  
limited guarantee

READ ALL INSTRUCTIONS BEFORE USE.  
SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

## IMPORTANT SAFEGUARDS :

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons. THIS IS NOT A TOY.
- DO NOT place or drop into water or any other liquid.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has been damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep away from heated surfaces.
- NEVER drop or insert any object into any opening.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, send to a HoMedics Service Centre. All servicing of this appliance must be performed by authorised HoMedics service personnel only.

## SAFETY CAUTIONS :


PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

- **If you have any concerns regarding your health, consult a doctor before using this appliance.**
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time of 15 minutes.

## RECOMMENDED CLEANING INSTRUCTIONS :

### Removable Cover

- Machine Wash Warm 40°, Do Not Iron, Tumble Dry Low.

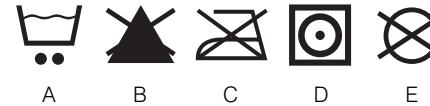- A Warm wash 40°
- B Do not bleach
- C Do not iron
- D Tumble dry at minimum setting
- E Do not dry clean

### Water Bottle

- Clean only with a soft, slightly damp sponge. Never allow water or any other liquids to come into contact with the battery compartment.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner to clean.
- Do not attempt to repair the Water Bottle. There are no user serviceable parts. For service, contact HoMedics Customer Service (see guarantee section for contact details).

## INSTRUCTIONS FOR USE :

### Removing the Massage Mechanism

1. Unscrew the massage mechanism using the outer lugs marked 'Lock - Unlock' (Fig 1 & 2).
2. Slide the massage mechanism out of the water bottle (Fig 3 & 4)

Fig 1

Fig 2

Fig 3

Fig 4

## Installing/Changing the Batteries

1. Unscrew the battery cap from the massage mechanism using the centre 2 lugs marked 'Lock – Unlock' (Fig 5 & 6).
2. Place 2 x AA batteries into the chamber (negative [-] first) as indicated on the compartment (observe +/- polarity as indicated on the outside of the massage unit and yellow label inside the compartment (Fig 7, 8 & 9).
3. Replace the battery cap on the massage mechanism using the centre 2 lugs marked 'Lock – Unlock' (Fig 10 & 11).

Fig 5

Fig 6

Fig 7

Fig 8

Fig 9

Fig 10

Fig 11

## Battery Precautions

- Use only the size and type of batteries specified (2 "AA").
- When installing batteries, observe proper +/- polarities. Incorrect installation of battery may cause damage to the unit.
- Do not mix different types of batteries together (e.g., alkaline with carbon-zinc or old batteries with new ones).
- If the unit is not going to be used for a long period of time, remove batteries to prevent damage due to possible battery leakage.
- Do not dispose of batteries in fire. Batteries may explode or leak.
- Do not immerse battery compartment in water.

## Filling the Water Bottle

1. Unscrew and remove the massage mechanism as described in 'removing the massage mechanism' section.
2. Slide the massage mechanism out of the water bottle (Fig 12).
3. Holding the water bottle open in an upright position (Fig 13).
4. When filling the hot water bottle, do not use boiling water and fill to a maximum of two-thirds capacity (1.25 ltrs, 42 fl. oz.) with hot water from a kettle which has been left to stand for one minute. Hold the bottle by the neck in an upright position and fill slowly to avoid hot water splashing back (Fig 14 & 15)
5. Insert the massage mechanism stopper to within 2cm of closing (Fig 16). Then expel the excess air by lowering it carefully onto a flat surface until water is just visible through the opening.
6. Screw the massage mechanism sufficiently tight to ensure that there is no leakage. Finger tight should be adequate (Fig 17).

Fig 12

Fig 13

Fig 14

Fig 15

Fig 16

Fig 17

## NOTE:

- Do not fill using water from the domestic hot water system as this can considerably shorten the life of the hot water bottle.
- Prevent contact with hot surfaces.
- Prevent contact with oil or grease.
- When not in use, drain completely and keep, with the stopper removed, in a cool, dry, dark place. Prevent exposure to sunlight. Do not place anything on top of the bottle.
- Check the bottle and stopper for wear and damage at regular intervals.
- Replace water bottle cap and make sure it is on securely
- Test for leakage by turning bottle upside down and pressing firmly.
- Water Bottle may be used without massage function.

**Caution:**

- Do Not overfill.
- Do Not use boiling water.

**Operating Instructions**

- To activate massage press blue power button (Figure 18).
- Hold Water Bottle over desired area to be massaged.
- After completing massage turn Water Bottle off by depressing the blue power button once.
- Empty water bottle when not in use, by unscrewing water bottle cap and removing it from water bottle. Pour water into the sink.

**Fig 18****Removing the Cover**

- To remove cover, first make sure water bottle is empty.
- Remove water bottle cap / massage tube.
- Pull hook and loop closure apart (Figure 19) and pull water bottle out through this opening (Figure 20).
- To replace cover, make sure water bottle is empty. Pull hook and loop closure apart and slide water bottle in bottom first (Figure 21). Bend top part and push the top through the top opening. Secure hook and loop closure together.

**Fig 19****Fig 20****Fig 21**

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