ELLE MACPHERSON THE BODY™

Cellulite Massager

2 Year Guarantee

Instruction Manual

By
HOMEDICS®

ELM-CELL100-EU

Thank you for purchasing the Cellulite Massager, the latest addition to my new 'The Body' range. My team and I have developed the range to help you get the very best results from your beauty regime. Intuitive and simple to use, it's never been easier to discover my professional beauty secrets.

No-one likes to talk about cellulite, but it can affect anyone, young and old, skinny or curvy. Finding an effective solution to uneven skin tone is not an easy task, but I know that cellulite massage really can help. My Cellulite Massager encourages lymphatic drainage and the removal of toxins, improving circulation and helping to break up the fatty deposits that can result in cellulite. I would recommend cellulite massage as part of a healthy and active lifestyle – a balanced diet, drinking plenty of water (I drink 3 litres a day), regular exercise and a targeted massage twice a week.

Elle 😊

**READ ALL INSTRUCTIONS BEFORE USE.
SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.**

IMPORTANT SAFEGUARDS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.

- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT use outdoors.
- DO NOT crush. Avoid sharp folds.
- This product requires a 220-240 V AC power.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. For service, send to a HoMedics Service Centre. All servicing of this appliance must be performed by authorised HoMedics service personnel only.
- Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch. Keep air openings free of lint, hair etc.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- Failure to follow the above may result in the risk of fire or injury.
- Burns can result from improper use.

SAFETY CAUTIONS :

PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

- **If you have any concerns regarding your health, consult a doctor before using this appliance.**
- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by Diabetics
- DO NOT use on an infant, invalid or on a sleeping or unconscious person. DO NOT use on insensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Do not use for longer than the recommended time.
- Care should be taken when using heated surfaces. If the product feels excessively hot, switch off at mains and contact the HoMedics Service Centre.
- You may soften massage force by placing a towel between yourself and the unit.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Please ensure that all hair, clothing and jewellery are kept clear of the massage mechanism or any other moving parts of the product at all times.

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

WHAT IS CELLULITE?


Cellulite is a condition in which the skin has a dimpled, lumpy appearance due to underlying fat deposits. It is most noticeable on the buttocks and thighs, but also occurs on the upper arms.

By using this product you can help to improve circulation, encourage lymphatic drainage, when combined with eating a healthful, balanced diet and exercising you can help your body to reduce the fat content in cells and reduce the appearance of cellulite.

PRODUCT FEATURES

Ergonomic and lightweight design

Comfortable grip
Fits the palm of your hand

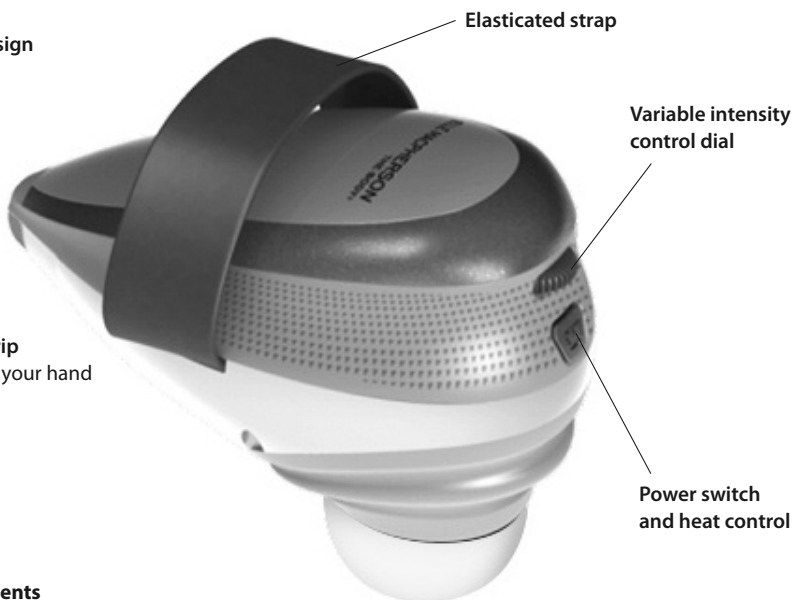Three Attachments

1. Dome – use for targeted massage relief
 2. Small cellulite attachment – Use for cellulite massage of smaller areas e.g arms
 3. Large cellulite attachment – Use for cellulite massage of larger areas e.g. thighs, buttocks
- Heat function can be activated with ALL attachments.

Fig. 1

INSTRUCTIONS FOR USE

Massage Function

NOTE: Always ensure the appliance is set to the 'OFF' before inserting the plug into a mains outlet.

1. Connect the power cable to the product by inserting the connector into the socket attached to the product. Plug the appliance into a 230V mains outlet and press the power button once to switch on.
2. Your massager comes equipped with 3 attachments (Figure 1). Attach the different heads by twisting on (clockwise) or twisting off (anti-clockwise) (Figure 2).
3. Grip massager by placing your hand through the elasticated strap and grasping the unit with fingers outstretched. The strap should be positioned over the back of your hand (Figure 3).
4. Press the power button once to activate the massager. The button will illuminate blue. Adjust the massage intensity by rotating the dial on the front of massager (Figure 4).
5. **Heat Function:** All attachments are compatible with the heat function. To activate press the power button a 2nd time, it will illuminate red to indicate that the heat setting on. The centre of the attachment will begin to heat up. **NOTE: To eliminate the risk of being burned, never use the heat function without the attachments. Do not use the thermal heat continuously for more than 15 minutes. Do not use the thermal heat if you have sensitive skin, poor circulation, discoloured areas, swelling, burning and inflammation or where skin eruptions or sores are present. Heat should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.**
6. Apply the massager head lightly to the area you wish to treat. Move the unit slowly, holding it over the targeted area. Treat the area for a short period of time and move to another area. NEVER treat the same area for longer than 3 minutes.
7. To turn the unit off, press the power switch a third time. For your safety the device will automatically turn off after 15 minutes of continuous use.

Tips for Use

Massage in a firm circular motion moving towards the heart, concentrating on targeted problem areas of the body.

The Elle Macpherson Cellulite Massager is especially effective, when used regularly and in conjunction with a healthy balanced diet and undertaking of regular exercise.

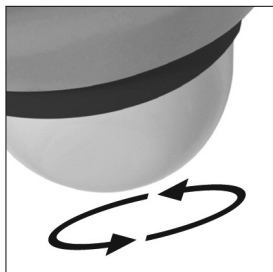
Fig. 2

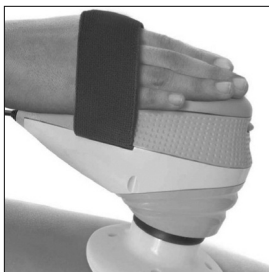
Fig. 3

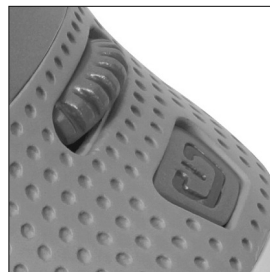
Fig. 4

MAINTENANCE:

Cleaning

Unplug the appliance and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.

- Never allow water or any other liquids to come into contact with the appliance
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

Storage

Place the appliance in its box or in a safe, dry, cool place. DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

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